


Universal Smartphone Connector for HEINE DELTA 30 and DELTAone. Digital dermatoscopic documentation using a smartphone camera.

IT'S JUST THAT EASY!



1 You can use almost any smartphone with the HEINE Universal Smartphone Connector, regardless of operating system (Android, Windows, ...). HEINE has phone cases in its product range that fit most iPhone* models, where the adaptor is already permanently attached to the HEINE phone case.

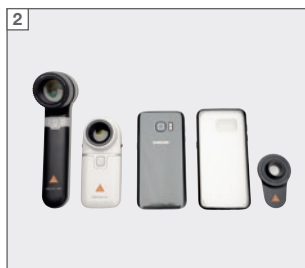

2 You need a DELTA 30 or DELTAone, the smartphone connector and a phone case that fits your smartphone**.




3 The connector is only suitable for phone cases with a flat surface; the camera should not protrude out of the case.



4 If your smartphone has just one camera, affix the connector directly over the centre of this camera.

If your smartphone has multiple cameras, first determine which of the cameras is the main camera: place the connector over the centre of one of the camera lenses on the back of the smartphone. Is the field of view on the screen clear in the camera view? If everything looks correct, affix the connector over this camera.

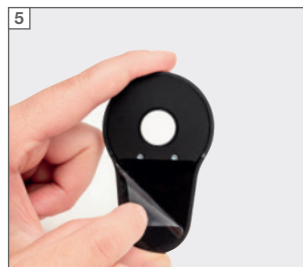

5 Remove the protective film and press the connector firmly over the centre of the selected camera lens.



6 Almost done!!



7 With the phone now in its case, the magnet on the connector keeps the phone in position on the focusing ring of the dermatoscope.

You can now use the camera function on your smartphone as normal to take photos.

*Apple, the Apple logo and iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. HEINE recommends using iPhones from version 6 upwards. All other trademarks are the property of their respective owners.

**The smartphone and phone case are not included in the scope of delivery.